

ASAHI

Peripheral Guide Wire

ASAHI's
composite guide wire




ASAHI

Peripheral Guide Wire

Exceptional guide wires for precise interventions.

One Piece Core Wire

Provides precise torque response and smooth device delivery.

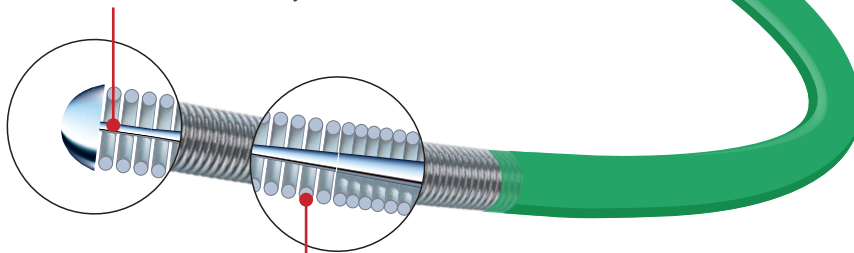


Joint-less Spring Coil

Ideal torque transmission with reliable trackability.

**Regalia
Treasure
Astato**



ACT ONE technology improves torque transmission and durability.

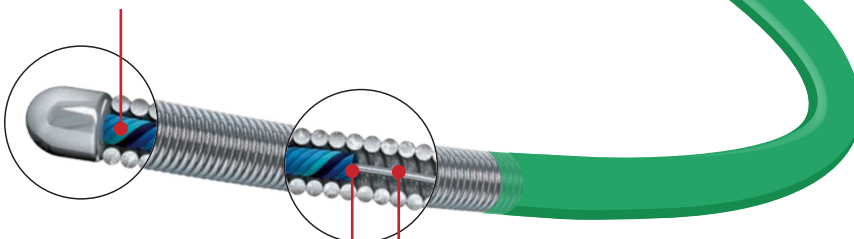

Round Core Wire

Improves torque response.

Balanced Support Shaft

Makes it easier to transmit push force from shaft to tip.


ASAHI **Gladius**
ASAHI **Halberd**
ASAHI **Gaia PV**

Frontline


Regalia XS 1.0 0.36mm (0.014inch) Sub-total/Tortuous




Balanced lubricity and trackability for BTK procedures.

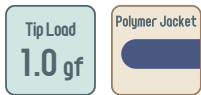

Treasure Floppy 0.45mm (0.018inch) Workhorse

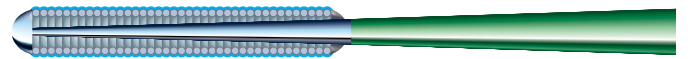

Precise handling with a flexible tip for tortuous anatomy.

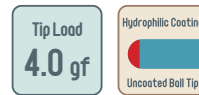

Specialty

ASAHI Halberd 0.36mm (0.014inch) Complex Lesion




Composite guide wire improves torque for navigation, and the Micro-cone tip improves penetration.




Treasure 12 0.45mm (0.018inch) Complex Lesion




Consistent response and tactile feedback for complex lesion navigation.

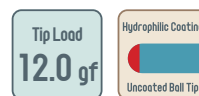

Astato XS 20 0.36mm (0.014inch) Highly Calcified Lesion

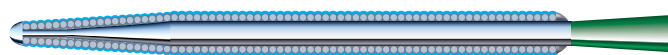


Tapered tip provides excellent penetration.

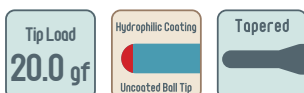

Astato 30 0.45mm (0.018inch) Highly Calcified Lesion



Tapered tip provides excellent penetration.

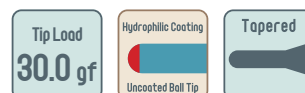

ASAHI Gladius 0.36mm (0.014inch) Workhorse




Support ● ● ○ ○ ○ Composite guide wire improves tip durability and torque to treat lesions in multiple vessels.

Torque ● ● ● ○ ○

Penetration ● ● ○ ○ ○



ASAHI Gladius 0.45mm (0.018inch) Workhorse

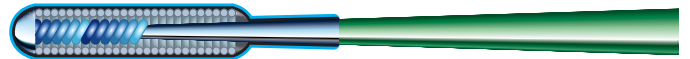


Support ● ● ● ● ● Composite guide wire improves torque response with increased support for device delivery.

Torque ● ● ● ○ ○

Penetration ● ● ○ ○ ○



ASAHI Halberd 0.45mm (0.018inch) Complex Lesion

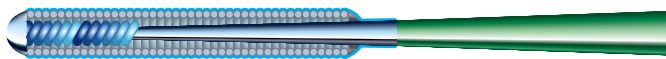


Support ● ● ● ● ○ Composite guide wire provides precise torque for navigation with more support and Micro-cone tip for proximal cap penetration.

Torque ● ● ● ● ●

Penetration ● ● ● ● ○



ASAHI Gaia PV 0.45mm (0.018inch) Complex Lesion

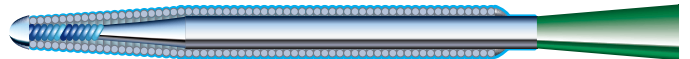


Support ● ● ● ● ○ Composite guide wire with a lower tip load improves navigation and torque. The tapered coil and Micro-cone tip makes it easier to advance.

Torque ● ● ● ● ●

Penetration ● ● ● ● ○




Astato XS 40 0.36mm (0.014inch) Highly Calcified Lesion

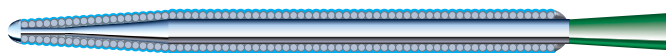


Support ● ● ○ ○ ○ Tapered tip provides excellent penetration.

Torque ● ● ● ● ○

Penetration ★ ★ ★ ★ ★

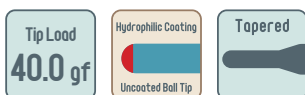

Wire Selection Chart

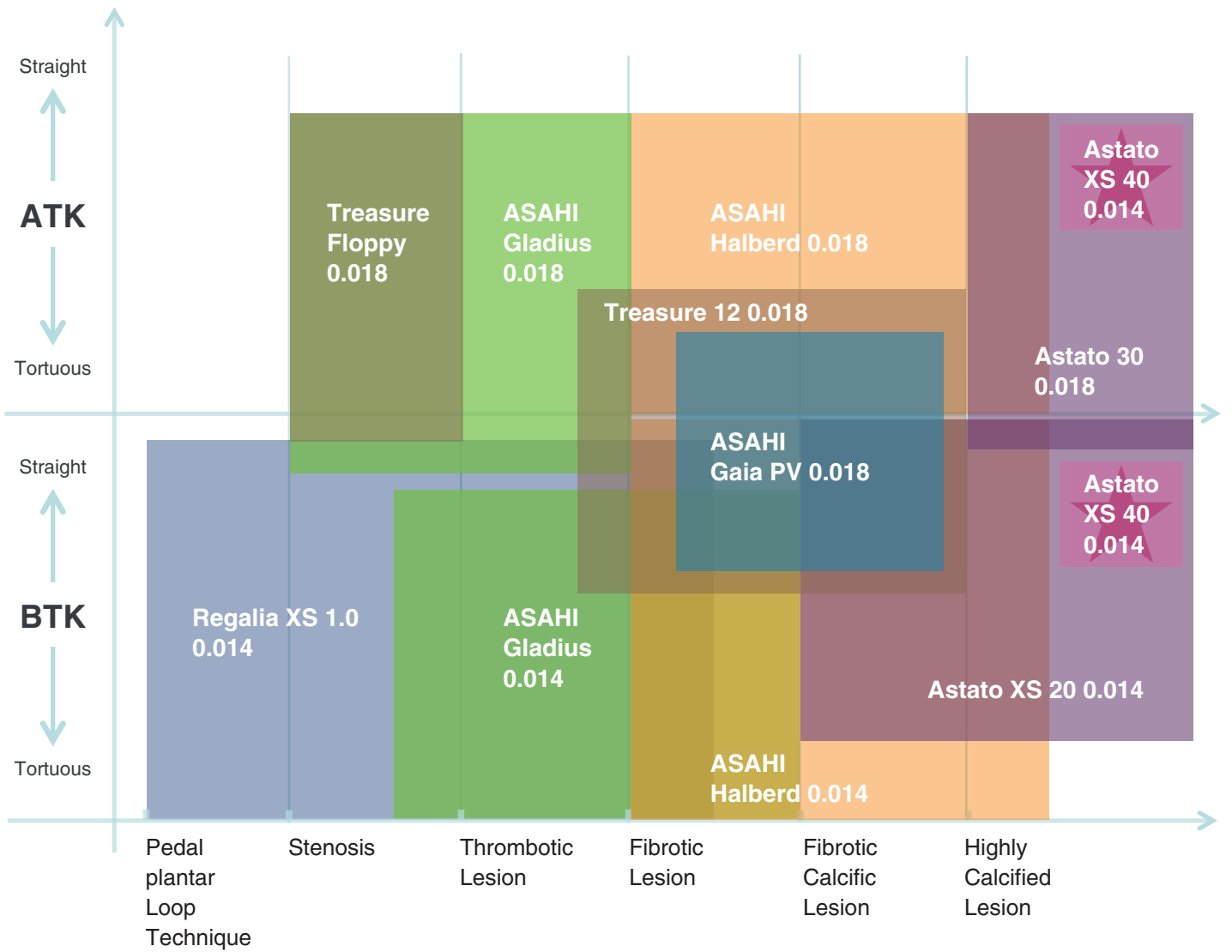

Ordering information

Product	Diameter	Catalog No.	Length	Tip Shape	Composite guide wire
ASAHI Gladius	0.36mm (0.014inch)	PPW14R100S	200 _{cm}	Straight	ACT ONE
		PPW14R200S	235 _{cm}		
		PPW14R300S	300 _{cm}		
	0.36mm (0.014inch)	PPW14R100P	200 _{cm}	Pre-shape	
		PPW14R200P	235 _{cm}		
		PPW14R300P	300 _{cm}		
	0.45mm (0.018inch)	PPW18R100S	200 _{cm}	Straight	
		PPW18R200S	235 _{cm}		
		PPW18R300S	300 _{cm}		
ASAHI Halberd	0.36mm (0.014inch)	PHW14R101S	200 _{cm}	Straight	ACT ONE
		PHW14R201S	235 _{cm}		
		PHW14R301S	300 _{cm}		
	0.36mm (0.014inch)	PHW14R101P	200 _{cm}	Pre-shape	
		PHW14R201P	235 _{cm}		
		PHW14R301P	300 _{cm}		
	0.45mm (0.018inch)	PHW18R101S	200 _{cm}	Straight	
		PHW18R201S	235 _{cm}		
		PHW18R301S	300 _{cm}		
	0.45mm (0.018inch)	PHW18R101P	200 _{cm}	Pre-shape	
		PHW18R201P	235 _{cm}		
		PHW18R301P	300 _{cm}		
ASAHI Gaia PV	0.45mm (0.018inch)	PHW18R102S	200 _{cm}	Straight	ACT ONE
		PHW18R202S	235 _{cm}		
		PHW18R302S	300 _{cm}		
	0.45mm (0.018inch)	PHW18R102P	200 _{cm}	Pre-shape	
		PHW18R202P	235 _{cm}		
		PHW18R302P	300 _{cm}		
Regalia XS 1.0	0.36mm (0.014inch)	PAGP140000	180 _{cm}	Straight	—
		PAGP140300	300 _{cm}		
Astato XS 20	0.36mm (0.014inch)	PAGH143092	180 _{cm}		
		PAGH143392	300 _{cm}		
Astato XS 40	0.36mm (0.014inch)	PAGHW143094	200 _{cm}		
		PAGHW143394	300 _{cm}		
Astato 30	0.45mm (0.018inch)	PAGH18M071	180 _{cm}		
		PAGH18M371	300 _{cm}		
Treasure Floppy	0.45mm (0.018inch)	PAGH18M072	190 _{cm}		
		PAGH18M372	300 _{cm}		
Treasure 12	0.45mm (0.018inch)	PAGH18M070	180 _{cm}		
		PAGH18M370	300 _{cm}		

Indications for Use

This Product is intended to facilitate the placement and exchange of diagnostic and therapeutic devices during intravascular procedures.
This device is intended for peripheral vascular use only.

Not all shapes and sizes are available in all markets. Please check availability with your sales representative.
Coated with SLIP-COAT® coating.



© 2018 ASAHI INTECC CO., LTD.
"ASAHI", "Astato", "ASAHI Gaia", "ASAHI Gladius" and "Treasure" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.

Manufactured by:

Distributed by:



ASAHI INTECC CO.,LTD.

3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan

Contact phone number: +81-561-86-9101

<http://www.asahi-intecc.com>